

Pull Cord Wall Button Quick Start Guide

When used in conjunction with the Cellular Base Station and/or Mobile Device, the Pull Cord Wall Button allows you to obtain emergency help with multiple methods; a press of a button, a pull of a cord, or by calling out for help*. The PCWB can be mounted on a flat surface, such as a wall, or placed on a table top. Mounting hardware and two-sided tape is included. The PCWB includes 2 C batteries that are easy to replace.

In an emergency, you may use any of the following forms to send a signal for help:

1. Press the Emergency Call Button
2. Pull the cord out
3. Shout out **"Call for HELP"** twice within 10 seconds



*Voice command may not trigger signal for help 100% of the time. If able, users should always push their Help Button when they need assistance. If installed in the shower, voice activation may be difficult due to background noise. Please stand within 10 feet of the device for best voice recognition.

MHP-06908A

Wall Installation

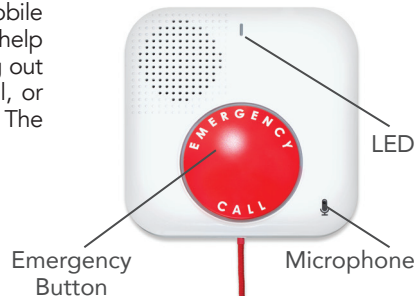


Table Installation



VAB-02

How to mount the PCWB on the wall:

- 1 Clean and dry the surface.
- 2 Peel the strips off the mounting tape and press the PCWB into place. (see **Figure A**)
- 3 Insert the pull cord until the red LED illuminates, indicating the cord has been fully inserted. (see **Figure B**)
- 4 The pull cord should only be used in a configuration without the table mount bracket.



Figure A



Figure B

How to install the PCWB on a table top:

Snap the PCWB into the table top bracket. (see **Figure C**)



Figure C



How to replace the batteries:

- 1 When the red LED flashes momentarily every 2 minutes, it is time to replace the batteries. Pull off the back cover. (see **Figure D**)
- 2 Insert the batteries and wait for the red LED to illuminate. (see **Figure E**)
- 3 Close the back cover. (see **Figure F**)



Figure D

Hold here and pull off the back cover



Figure E

Align white post with black hole



Figure F

VAB-02