

Mobile Medical Alert Systems

Complete Protection, In-Home Wireless and Mobile On-the-Go Systems
(Optional Internal Fall Detection on Mobile Device)



Complete User Guide

IMPORTANT NOTICE

Internal Fall Detection does not detect 100% of falls. If able, users should always push the Help Button when they need assistance. Mobile Device is not intended to replace a caregiver.

Fall Detection: If you added our optional Fall Detection service during sign up, please note that detecting falls works when the Mobile Device is worn around your neck. Failure to wear the Mobile Device with Fall Detection around your neck will result in missed falls and misinterpretation of falls. Wear the Mobile Device on the outside of your clothing as to not muffle the microphone and speaker allowing you to communicate with the central station operator more clearly. Additionally you will have better access to press the emergency button when you need assistance. Please remember to handle Mobile Device with care when putting it on and taking it off in order to avoid false alarms.



If you do not have the optional Fall Detection service on your Mobile Device and would like to add it, please call our Customer Support team using the phone number listed on the contact information insert that came with your system to find out how you can add Fall Detection service for a low price in addition to your system plan!

Note: Please handle the Mobile Device with care when putting it on or taking it off, as it may interpret this movement as a fall and activate. If it does, the red LED will begin flashing and you will hear the announcement, "Fall detected, press and hold button to cancel." You can cancel the alarm by pressing and holding the emergency button approximately 5 seconds, until you hear the announcement "Alarm Cancelled." If you are unable to cancel, please tell the Emergency Operator that "THIS IS A FALSE ALARM." If you do not explain that this is a false alarm, EMS will be contacted.

If you have any questions, please call Technical Support.

SRDR-00025 Rev. 4

Complete User Guide

This User Guide includes information covering the 4 most popular systems:



**Complete Protection System
(DUO)**



**In-Home Wireless System
(Classic)**



**Mobile On-the-Go System
(SOLO)**




**Mobile Double System
(Mobile DUO)**

Note: Remember to always carry the Mobile Device with you when going outside of your home.

SRDR-00025 Rev. 4

Complete Protection System

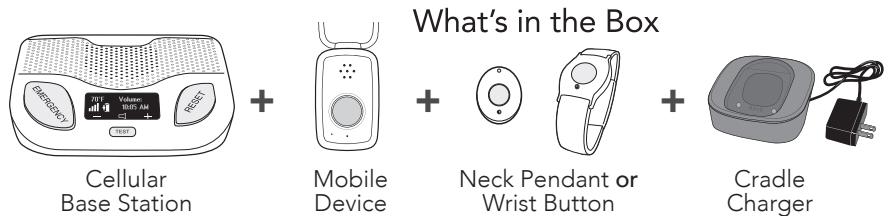


Optional Accessories




The Fall Button™, Wall Button and Mobile Case are optional and at an additional cost. The Fall Button™ has a separate user guide. **Please note:** there is a difference between optional internal fall detection on your mobile device and a Fall Button™.

SRDR-00025 Rev. 4

In-Home Wireless System

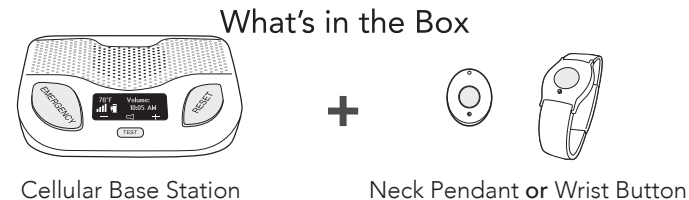

Optional Accessories

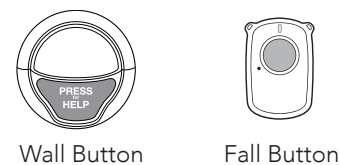


The Fall Button™, and Wall Button are optional and at an additional cost. The Fall Button™ has a separate user guide. **Please note:** there is a difference between optional internal fall detection on your mobile device and a Fall Button™.

SRDR-00025 Rev. 4

Mobile On-the-Go System

What's in the Box

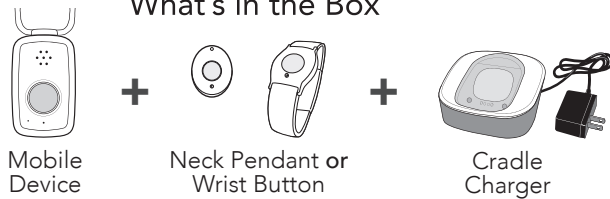

Optional Accessories




The Fall Button™, Wall Button and Mobile Case are optional and at an additional cost. The Fall Button™ has a separate user guide. **Please note:** there is a difference between optional internal fall detection on your mobile device and a Fall Button™.

SRDR-00025 Rev. 4

Mobile Double System

What's in the Box

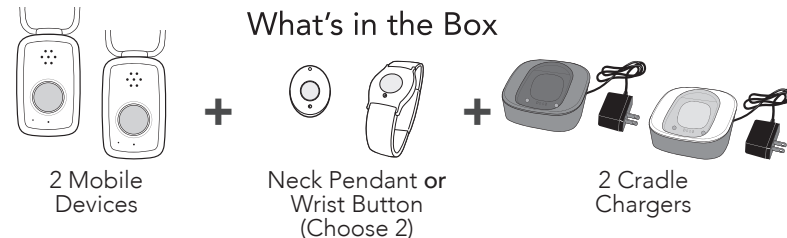

Optional Accessories



The Fall Button™, Wall Button and Mobile Case are optional and at an additional cost. The Fall Button™ has a separate user guide. **Please note:** there is a difference between optional internal fall detection on your mobile device and a Fall Button™.

SRDR-00025 Rev. 4

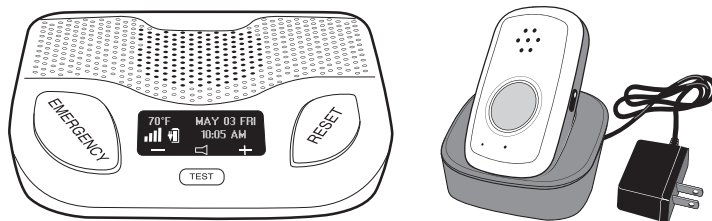


Figure 1

Setup Guide for the Medical Alert System

We thank you for choosing a Medical Alert System and welcome you to the family.
Congratulations on making a smart choice!

If you have any questions during the setup process, please call Technical Support.

SRDR-00025 Rev. 4

Cellular Base Station Features (Front)

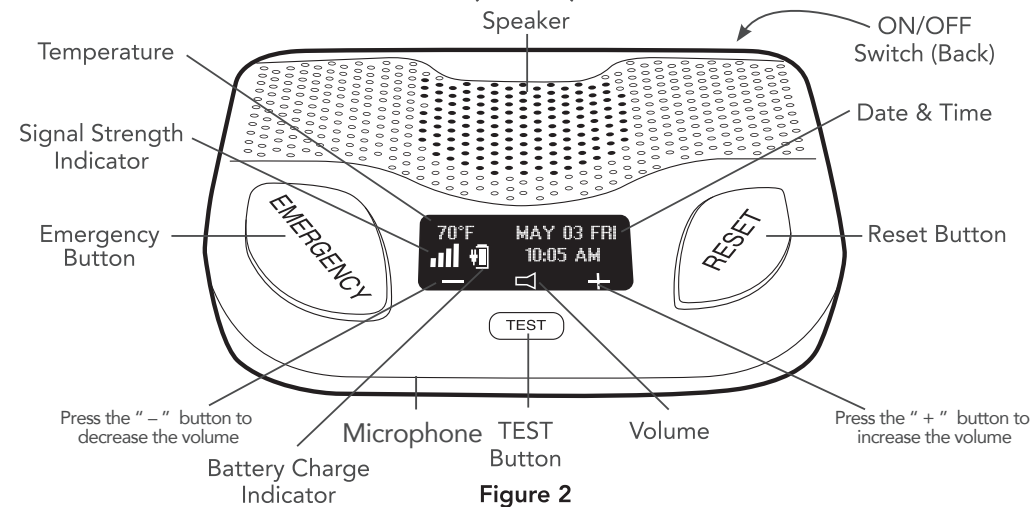


Figure 2

SRDR-00025 Rev. 4

Mobile Device System Features


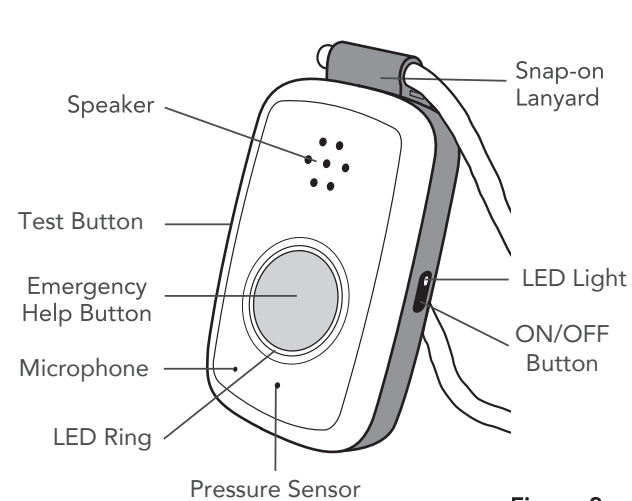


Figure 3



SRDR-00025 Rev. 4

Waterproof Help Buttons

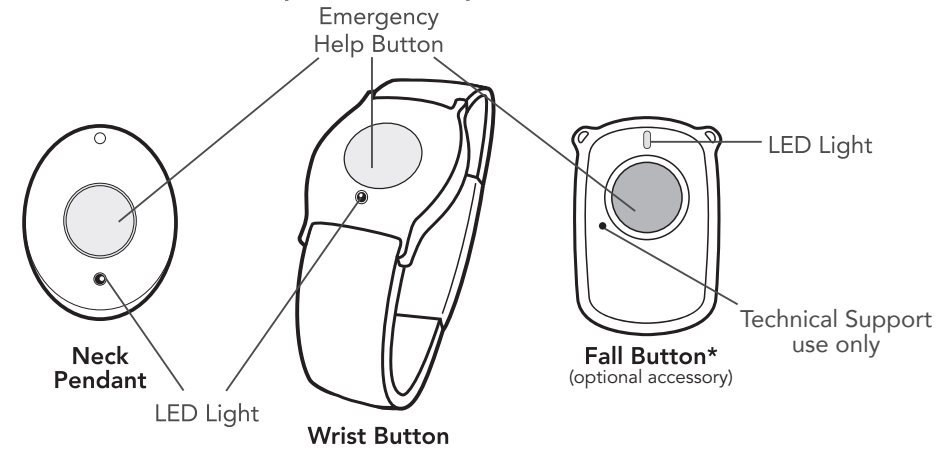

*Fall detection does not detect 100% of falls. If able, users should always push their Help Button when they need assistance. Fall Button is not intended to replace a caregiver.

SRDR-00025 Rev. 4

Please follow these simple steps to get your system set up and tested.

Connecting your Base Station

Step 1. Select a Location

- We recommend you place the Base Station near the center of your home, such as your living room, family room or kitchen area.
- Place it near an electrical outlet.
- Do not place any object in front of the Base Station.
- Do not place near any appliances that make noise.

Step 2. Plug in the Power Cord

(Please refer to figure 4)

Plug the power cord that is attached to the back of the Base Station into an electrical outlet.

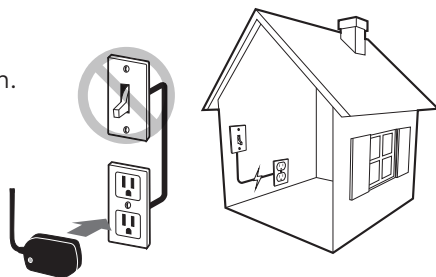


Figure 4

To avoid accidentally turning off the Base Station, **DO NOT** plug it into an electrical outlet that is controlled by a light switch.

SRDR-00025 Rev. 4

Step 3. Turning ON your Base Station

- Turn on your Base Station using the ON/OFF switch located on the back. (Please refer to figure 5)
- The display screen will turn on. (Please refer to figure 6)
- The Emergency and Reset Buttons will illuminate.
- If the Base is setup correctly, the Base Station will say "system ready" within 60 seconds.
- Make sure your Base Station is connected to the cellular network by checking for the signal strength bars on the display screen.
- You can use the volume control buttons on the front of the unit at the bottom of the display screen, indicated by the "-" or the "+" sign, to adjust if the voice is too loud or too soft.
- Within a few minutes of turning on the system, the Cellular Base Station will obtain the date and time from the cellular network and they will appear in the display area. You are not able to set the date and time manually.

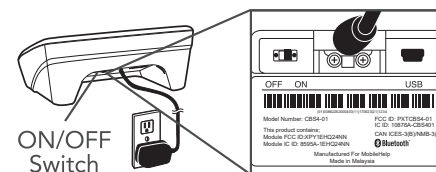


Figure 5

SRDR-00025 Rev. 4


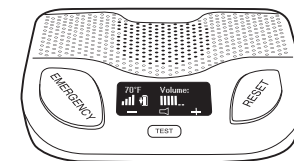


Figure 6

Setting up your Mobile Device

Step 1. Select a Location

- We recommend you place your Mobile Device in your bedroom or any other place around your home where it can be easily accessible.
- Place your Mobile Device's cradle charger near an electrical outlet.

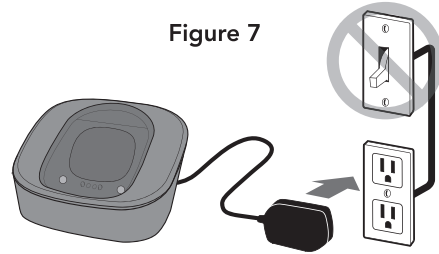

Step 2. Connect the Cradle Charger (Please refer to figure 7)

- Plug the AC adaptor into an electrical outlet that is not controlled by a light switch.



WARNING

To prevent electrical shock please keep the system away from wet locations.

AVERTISSEMENT

Pour éviter le choc électrique, SVP tenir le système éloigné des endroits humides.

SRDR-00025 Rev. 4

Step 3. Charge your Mobile Device

Charge your Mobile Device every night to ensure a fully charged battery during waking hours. Placing the Cradle Charger near the bed will make it easier to speak to an operator if you need to press the emergency button, or pendant, during the night.

- Place the Mobile Device in the cradle charger. (Please refer to figure 8)
- The LED light on your Mobile Device may display multiple colors when first placing it on the charging cradle.
- The Mobile Device will turn on within a few seconds and the LED indicator around the emergency button will illuminate. It normally takes 30-60 seconds for the Mobile Device to register with the cellular network.
- When the Mobile Device is setup correctly and connects to the cellular network, it will announce "system ready."
- The LED indicator around the emergency button will illuminate red, indicating that the device still needs charging. When the device is fully charged, the LED indicator will change from red to green and will remain green until the device is removed from the charging cradle. It will take up to 3 hours for the device to be fully charged.

SRDR-00025 Rev. 4

To turn on your Mobile Device:

- **Option 1:** Simply place the Mobile Device in the cradle charger, assuming the cradle charger is plugged in.
- **Option 2:** Press the ON/OFF button on the side of the Mobile Device and hold for approximately 3 seconds. (Please refer to Figure 8)
- All the lights will turn on and the green LED light on the power button will start blinking. This indicates the Mobile Device was successfully turned on and it will announce "system ready."

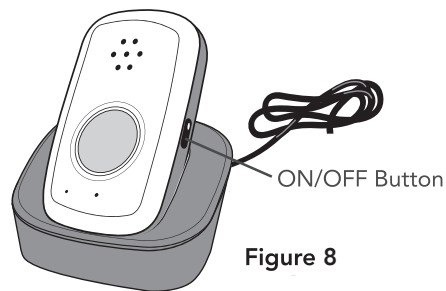


Figure 8

To turn off your Mobile Device:

- Press and hold the power button for about 3 seconds. The device will emit a beep and the power button LED no longer blinks and all LEDs remain off.

SRDR-00025 Rev. 4

Testing your Complete Protection System

Please test your system on a weekly basis using the TEST button.

- Always have your Base Station and Mobile Device nearby when testing your system. It is possible that the alarm may go to a live operator.
- You can use the TEST button on your Base Station and Mobile Device to conduct a test without speaking to an Emergency Response Operator by following these steps:
 - Press and hold the TEST button on Base Station until it illuminates green.
 - Press and hold the TEST Button on Mobile Device (the button on the side with the letter "T").
- You will hear an announcement over the Base Station and Mobile Device that will instruct you to press the Emergency Button or Pendant.
- Press the Neck Pendant or Wrist Button or Emergency Button on the base or Mobile device.
- The Base Station and Mobile Device will announce "test call sent to Emergency Response Center."

SRDR-00025 Rev. 4

- If your test was successful, you will hear “thank you for testing your device” from the Base Station and Mobile Device.
- If your test was not successful, the Base Station or Mobile Device will announce “user auto-test failed.” Please contact Technical Support.
- If by chance the operator speaks to you, please state clearly that “THIS IS A TEST.” If you do not explain to the operator that you are okay and that you are testing, EMS will be contacted.
- **The Base Station contains a battery backup. If you would like to test the battery backup capability, please unplug the power for the Cellular Base Station from the AC outlet. Then follow the testing steps above. Once you have finished testing, please remember to plug the power back into the AC outlet.**

SRDR-00025 Rev. 4

Testing your In-Home Wireless System

Please test your system on a weekly basis using the TEST button

You can use the TEST button on your Cellular Base Station to conduct a test without speaking to an Emergency Response Operator by following these steps:

- Press and hold the TEST button on Base Station until it illuminates green.
- You will hear an announcement over the Base Station that will instruct you to press the Emergency Button or Pendant.
- Press the Neck Pendant or Wrist Button.
- The Base Station will announce “test call sent to Emergency Response Center”.
- If your test was successful, you will hear “thank you for testing your device” from the Base Station speaker.
- If your test was not successful, the Base Station will announce “user auto-test failed.” Please contact Technical Support.
- **The Base Station contains a battery backup. If you would like to test the battery backup capability, please unplug the power for the Cellular Base Station from the AC outlet. Then follow the testing steps above. Once you have finished testing, please remember to plug the power back into the AC outlet.**

SRDR-00025 Rev. 4

Testing your Mobile On-the-GO or Complete Protection On-the-GO Systems

For Mobile Device Light Indicator information, please see page 39.

MOBILE DUO CUSTOMERS:

If you have a Mobile DUO system, it will include 2 Mobile Devices. Please use the TEST button on both devices as detailed below.

- Press and hold the TEST button on Mobile Device (the button on the side with the letter "T"). (Please refer to figure 9)
 - You will hear an announcement over the Mobile Device that will instruct you to press the Emergency Button or Pendant.
 - Press the Neck Pendant, Wrist Button OR press and hold the Emergency Button on the Mobile Device for 2 seconds.
 - The Mobile Device will announce "test call sent to Emergency Response Center."
 - If your test was successful, you will hear "thank you for testing your device" from the Mobile Device.
- If the help button on your Mobile Device is not working properly or if the light is blinking rapidly, please call Technical Support.

Please remember, it is important to test your system weekly.

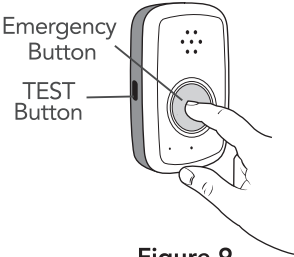


Figure 9

SRDR-00025 Rev. 4

Base Station Voice Guidance

Your Base Station has voice announcements. Select announcements are listed below:

Event	Voice Prompt	Notes
The Cellular Base Station is powered on	"System Ready"	Base Station ON/OFF switch has been turned on, be sure to connect to electric outlet or it will run on battery backup.
Both red Emergency and blue Reset Buttons are flashing	"Power not detected" or "Power detected" (when electricity is reconnected)	Check that the Base Station is correctly plugged into an electrical outlet.
The battery needs charging	"Low Battery"	Check that the Base Station is correctly plugged into an electrical outlet or it will run on battery backup.
The Emergency call is initiated	"Call in progress"	Alarm is being sent to the Emergency Response Center.
Call connected	"Please stand by for operator"	Respond to the emergency operator's call over the device speaker.
Cellular Base Station receives a Fall Signal	"Fall Detected - Press and hold Button to Cancel"	This will occur for about 20 seconds and will allow you time to cancel the alarm. During this time, you may press and hold the button on the Fall Button for 5 seconds to cancel the signal.
Fall Alarm is canceled	"Alarm Canceled"	Alarm has been canceled.

SRDR-00025 Rev. 4

Medication Reminders (optional service):

- At the scheduled time for your Medication Reminder, your Cellular Base Station will announce "Medication Reminder, please press reset to accept" and the blue Reset Button will flash.
- The name of the medication will be displayed on the Cellular Base Station screen, just as it was entered in the Connect web portal.
- The Base Station will make the announcement once a minute for 30 minutes or until the Reset Button is pressed.
- If you press the Reset Button, the reminder will be recorded as accepted. If you do not press the Reset Button within 30 minutes, the reminder will be recorded as ignored.
- For complete instructions for Medication Reminders, please log in to Connect web portal and under the "Help" tab, click on "Connect User Guide."
- If you have any question, please call Technical Support.

SRDR-00025 Rev. 4

Mobile Device 5.0

Voice Guidance

The Mobile Device 5.0 provides voice announcements. Select announcements are listed below:


Event	Voice Prompt
Mobile Device is powered on	"System Ready"
The emergency call is initiated	"Call in progress"
The battery is low	"Your battery is low. Please connect to a power source"
Mobile Device receives a Fall Signal	"Fall Detected - Press and Hold Button to Cancel"
Fall Alarm is canceled	"Alarm Canceled"
20-second window for cancellation has ended	"Fall Detected - Contacting Emergency Response Center"

Note: There may be situations where the voice prompts don't complete the entire phrase because another event occurs. In this instance, the latter voice prompt event will cut in and the current voice prompt may be terminated prematurely.

Radio Frequency

All devices frequency transmission = 433 MHz

SRDR-00025 Rev. 4

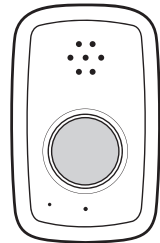

Fall Detection with Mobile Device 5.0

Voice Guidance

If your system includes the Fall Button or Internal Fall Detection, your Mobile Device will have these voice announcements listed below:



Event	Voice Prompt	Notes
Mobile Device receives a fall signal	"Fall detected. Please press and hold button to cancel"	This will occur for about 20 seconds and will allow you time to cancel the alarm. During this time, you may press and hold the button on the Fall Button for 5 seconds to cancel the signal.
Fall alarm is canceled	"Alarm canceled"	Alarm has been cancelled
20-second window for cancellation has ended	"Fall detected. Contacting emergency response center"	The signal will then be sent to the Monitoring Center and an operator will call your Mobile Device and/or the Base Station to assess your situation and dispatch the appropriate help to you.

Radio Frequency

All devices frequency transmission = 433 MHz

SRDR-00025 Rev. 4

Important Safety Information For All Systems

- You must always carry the Mobile Device with you when going outside of your home.
- Keep electrical cords away from heat sources and sharp edges.
- Keep electrical cords away from areas where someone could trip over them.
- Use caution when handling electrical cords.
- Never install electrical cords during a lightning storm.
- Do not use any attachment or accessory that is not intended for use with this system.
- Test your system using the TEST button when possible, as instructed in this User Guide.
- Do not block the Base Station's microphone. This will affect your ability to communicate with the Emergency Response Center.
- During a power outage, the Base Station's battery backup will typically last 24 hours after it has been charged for 3 hours.
- The Wrist Button and/or Neck Pendant are designed to work at an approximate range of up to 350 feet from the Mobile Device or up to 1400 feet from the Base Station, depending on the size and construction of your home.
- The Mobile Device can be used in the shower. Do not submerge in water.

Service availability and access/coverage on the cellular network is not available everywhere and at all times. Current GPS location may not always be available in every situation.

SRDR-00025 Rev. 4

- The Mobile Device contains sensitive hardware, take care when handling and avoid dropping it.
- Do not expose the Mobile Device to dust, dirt, sand or food.
- The battery is designed to last up to 3 days after it has been charged for 3 hours. Low cellular signal strength, no cellular signal, multiple button presses and optional internal fall detection service are factors which reduce battery life and necessitate charging of the mobile device sooner.
- The Mobile Device has a GPS receiver that works best when it has an unobstructed view to the sky.
- Wear your waterproof Help Button at all times.
- The Neck Pendant, Wrist Button and Fall Button are waterproof and can be taken with you in the shower. The Mobile Device is water resistant and can also be taken with you in the shower.
- Do not put the Neck Pendant or Wrist Button through the clothes washer or dryer. Do not attempt to dry your button out in the microwave or oven.
- If the LED on your Neck Pendant or Wrist Button flashes amber, the battery is low. Please contact MobileHelp Technical Support immediately to provide you with a replacement pendant.

WARNING

Strangulation and choking hazard. The lanyard has been designed to break away when tugged. However, the user can still suffer serious personal injury or death if the cord becomes entangled or stuck on objects.

AVERTISSEMENT

Risque de suffocation ou d'etranglement. La laniere est conçue pour se detacher lorsqu'elle est tiree, cependant, l'utilisateur peut subir une blessure ou la mort si la laniere se mele ou s'accroche a d'autres objets.

SRDR-00025 Rev. 4

Fall Detect Button and Internal Fall Detection FAQ (Optional Feature)

What happens if I fall?

If you fall, a signal will automatically be sent for help. If you are able to, you should always press the Emergency Button if you need help, as the Internal Fall Detection does not detect 100% of falls.

DETAILED PROCESS: Within a few seconds after a fall, the Internal Fall Detection will enter a 20 second CAUTION phase. The Mobile Device will start the announcement "Fall detected, press and hold button to cancel" and repeat it. The CAUTION phase is an early warning to cancel before help is contacted. After the 20 seconds, the Mobile Device enters the ALARM phase and announces, "fall detected, contacting emergency response center."

How can I cancel a Fall Detect alarm?

Alarms can be manually cancelled by pressing and holding the Emergency Button for at least 5 seconds during the CAUTION phase described on page 18.

If the alarm is not cancelled, an Emergency Operator will communicate with you through the Mobile Device. Please tell the operator that "THIS IS JUST A TEST." If the operator cannot hear you or you don't respond, EMS will be contacted.

If I fall and stand up, will the button still call for help?

If you fall, but then proceed to start moving around again, the button may still call for help. Many variables are considered in fall detection, for example how hard the fall is or how quickly you get back up.

SRDR-00025 Rev. 4

Frequently Asked Questions

How do I call for help?

Press your Neck Pendant, Wrist Button, or red Emergency Button on your Base Station or press and hold the emergency button on your Mobile Device for 2 seconds.

Do I need my Mobile Device when I leave my home?

Yes. Always take your Mobile Device with you when leaving your home. Wear your Neck Pendant or Wrist Button at all times, both at home and away, so help is never out of reach.

What happens if I push the button but I can't speak?

If you are unable to speak or be heard as a result of a fall, we will assume it is an emergency and notify Emergency Response on your behalf. We will also notify the personal contacts listed on your profile.

If I accidentally set off the alarm, what do I do?

If you accidentally set off the alarm, simply notify the Emergency Operator when they respond that "THIS IS A FALSE ALARM." The operator will disconnect and no further action will be taken.

How will an Emergency Response team get into my home?

Rescue personnel follow different procedures. In most cases, emergency personnel are authorized to enter your home to help you and cannot be held responsible for any damages caused as a result of gaining entry.

SRDR-00025 Rev. 4

We recommend that you purchase a "lock box" to safely store your house key. Your lock box code will be stored on your account and given to emergency personnel in order to gain entry to your home without causing damage. If you have your own lock box, please call Customer Service to ensure your code is noted on your account.

How far away from the Base Station will my Help Buttons work?

The Wrist Button and/or Neck Pendant are designed to work at an approximate range of up to 350 feet from the Mobile Device or up to 1400 feet from the Base Station, depending on the size and construction of your home.

What happens if I push my Help Button but cannot get to the Base Station?

If you push your Help Button, stay where you are. An Emergency Operator will attempt to communicate with you through your Base Station or Mobile Device. If you are not heard, we will assume it is an emergency and notify Emergency Response on your behalf.

Why is the Emergency Operator having trouble understanding me?

The microphones in the Base Station and Mobile Device are very sensitive, and will pick up surrounding noises from your television, radio, wind from ceiling fans, and the like. To increase the quality of the call, look for a quieter area to set up your Base Station or go to a quieter area with your Mobile Device.

SRDR-00025 Rev. 4

Why is my Help Button not activating my Base Station/Mobile Device?

Be sure to press the Help Button once firmly in the center to activate your Base Station and/or Mobile Device. A small red light below the Neck Pendant or Wrist Button will flash. If the Base Station and/or Mobile Device does not respond, check to ensure that the Help Button is in range of the Base Station and/or Mobile Device. If you still have issues with your system, please call Technical Support.

Is the Mobile Device waterproof?

The Mobile Device is water resistant and can be taken with you in the shower. The Help Buttons (Neck Pendant and Wrist Button) are waterproof and can be submersed in water.

What if I move to a different home, or my personal information has changed?

Please contact Customer Service any time there are any changes to your personal information, including changes to your address and phone number. Just give us a call on week before you move to determine whether the service will be available in your area.

What if I spend part of the year in another home?

We can transfer your service to your new location. Just give us a call one week before you go to determine whether the service will be available at your other home.

SRDR-00025 Rev. 4

I am going to travel. How should I prepare?

Your Mobile Device will work nationwide anywhere there is cellular network coverage. We recommend you call Customer Service at least 2 to 3 days prior to leaving, so we can update our systems with your new location address, emergency contacts, and lock box location. This ensures that if we have to alert emergency responders, they will have all relevant information to assist you without delays. Please remember to bring the following items: Neck Pendant and/or Wrist Button, Mobile Device, and Cradle Charger.

How do I prepare my device for air travel?

Your Mobile Device will need to be turned OFF. Simply press and hold the power button for about 3 seconds. The device will emit a beep and the power button LED no longer blinks and all LEDs remain off.

Is the Mobile Device lanyard adjustable?

The lanyard is non-adjustable. If you require an adjustable lanyard please call customer service to place your order.

SRDR-00025 Rev. 4

Can the Mobile Device cause interference with a pacemaker or other medical devices?

If you have a pacemaker or any medical device, consult your manufacturer for information on any possible risk of interference with the device. Our standard help buttons do not interfere with pacemakers.

Can I replace the Neck Pendant cord?

Yes, the Neck Pendant will work with just about any chain or cord, so you are able to use any one of your personal chains or necklaces. However, the risk of choking may increase if you do not use the lanyard provided.

Can I replace the Wrist Button band?

Yes, you may replace the Wrist Button band with a standard watch band of your choice.

Can I speak into my Help Button?

No, you can communicate with the monitoring center through your Base Station and/or Mobile Device only. Your Help Buttons **do not** have a speaker or microphone.

Does the Wrist Button strap contain any latex?

No, the strap is silicon based with no latex content.

FCC Statement

Product	Base	Neck Pendant	Wrist Button	Fall Detect Pendant	Wall Button
FCC ID	PXTCBS4-01	PXTIGP-01	PXTIGP-02	PXTPFD-01, PXT-RFDP-01	PXTWMP-01
FCC ID	PXTCBS4-02	PXTIGP-01	PXTIGP-02	PXTPFD-01, PXT-RFDP-01	PXTWMP-01
FCC ID	PXT-DC-MD5CCR-01	PXTIGP-01	PXTIGP-02	PXTPFD-01, PXT-RFDP-01	PXTWMP-01
FCC ID	PXT-DC-MD5-01	PXTIGP-01	PXTIGP-02	PXTPFD-01, PXT-RFDP-01	PXTWMP-01

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and (2) This device must accept any interference received, including interference that may cause undesired operation.

15.105(b)

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

15.21

You are cautioned that changes or modifications not expressly approved by the party responsible for compliance could void the user’s authority to operate the equipment.

RF Exposure Information

This device meets the government’s requirements for exposure to radio waves. This device is designed and manufactured not to exceed the emission limits for exposure to radio frequency (RF) Energy set by the Federal Communication Commission of the U.S. Government. This device complies with FCC radiation exposure limits set forth for an uncontrolled environment.

Cellular Base Station 4.0
Bluetooth Modem 2400 Mhz to 2483 Mhz

This equipment complies with IC RSS-102 RF exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with minimum distance 20cm between the radiator and your body.

.....
Déclaration d’exposition aux radiations:
Cet équipement est conforme aux limites d’exposition aux rayonnements IC établies pour un environnement non contrôlé. Cet équipement doit être installé et utilisé avec un minimum de 20 cm de distance entre la source de rayonnement et votre corps.

- 1. Canadian Compliance Statement**
This device complies with Industry Canada license-exempt RSSs. Operation is subject to the following two conditions:
1) This device may not cause interference, and
2) This device must accept any interference, including interference that may cause undesired operation of the device.

- Le présent appareil est conforme aux CNR d’Industrie Canada applicables aux appareils radio exempts de licence. L’exploitation est autorisée aux deux conditions suivantes :
- 1) l’appareil ne doit pas produire de brouillage;
 - 2) l'utilisateur de l'appareil doit accepter tout brouillage radioélectrique subi, même si le brouillage est susceptible d'en compromettre le fonctionnement.
2. This product meets the applicable Innovation, Science and Economic Development Canada technical specifications.
Ce produit répond aux spécifications techniques applicables à l’innovation, Science et Développement économique Canada.

Technical Information

Please refer to the website listed on the contact information card that came with your system.

Radio Frequency (RF)

All devices frequency transmission = 433MHZ

Notes

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

.....

Notes

This image shows a full page of white paper with horizontal dashed lines, typical of primary school writing paper. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

SRDR-00025 Rev. 4

Mobile Device Light Indicator

Device Status	LED Ring on Emergency Button
Power on	All lights (green, yellow, red) will turn on for 20-24 seconds
Pressing emergency button	Red LED ring will illuminate and will stay on until call has ended or after 2 minutes – whichever comes first. If in low battery status, red LED ring will continue blinking.
Low battery alert	Blinking red LED ring - you have about one hour of battery left
Charging	Solid red LED ring
Fully Charged	Red LED ring turns off and the green LED ring turns on, until the device is removed from the charge cradle
No cellular network	Blinking green LED ring
No cellular network while charging	Top half of the LED ring remains red to indicate charging while the bottom half of the LED ring blinks green to indicate no cellular network

SRDR-00025 Rev. 4

PLEASE REFER TO THE TECHNICAL SUPPORT NUMBER ON THE INSERT
THAT CAME WITH THIS GUIDE FOR ANY QUESTIONS OR CONCERNS.